

Challenger Ti-P

A global innovation in clip applier technology



Aesculap Surgical Technologies — Endoscopy

AESCULAP®

Challenger Ti-P

GLOBAL INNOVATION

The world's first **pneumatic** driven multi-fire clip applicator


- Easy and simple handling
- Fast, precise and reliable clip application
- Superior quality ligation clip
- Reusable applicator and single-use clip cartridge
- Handle is compatible with 5 mm and 10 mm shafts



Multi-Fire Ligation Clips


Challenger Ti-P multi-fire ligation clips are available in **Medium Large** and **Small Medium** sizes.

Each sterile package contains 1 clip cartridge and 1 CO₂ cylinder. The last clip in the cartridge is color-coded blue.

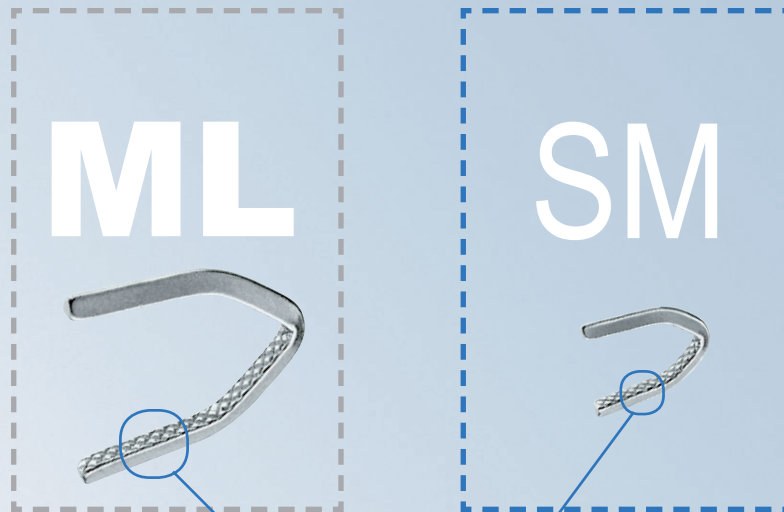


The unique inner clip texture maximizes clip adherence to vessels and helps prevent slippage.

ML
medium large

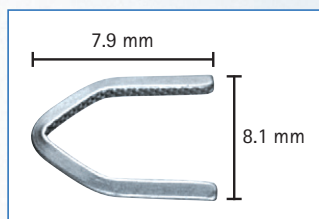

SM
small medium

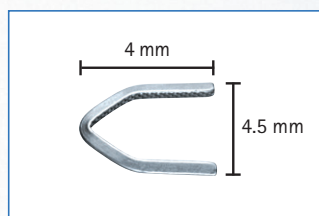


- Special closing mechanism allows clip to be repositioned prior to complete closure.



- When fully closed, the clip helps prevent tissue slippage.

Challenger Ti-P



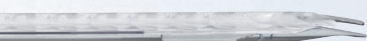

ML

Clip Applier

PL606R

Diameter = 10 mm
Working length = 370 mm

PL608R

Diameter = 10 mm
Working length = 260 mm

Clip Box



PL579T



- 12 sterile clip cartridges per box
-8 titanium clips per cartridge
-1 CO₂ cylinder per cartridge
- **Medium-large (ML) Clips**
for use with 10 mm clip appliers

SM

Clip Applier

PL604R

Diameter= 5 mm
Working length = 310 mm

Clip Box




PL574T



- 12 sterile clip cartridges per box
-12 titanium clips per cartridge
-1 CO₂ cylinder per cartridge
- **Small-medium (SM) Clips**
for use with 5 mm clip applier

Components and Spare Parts

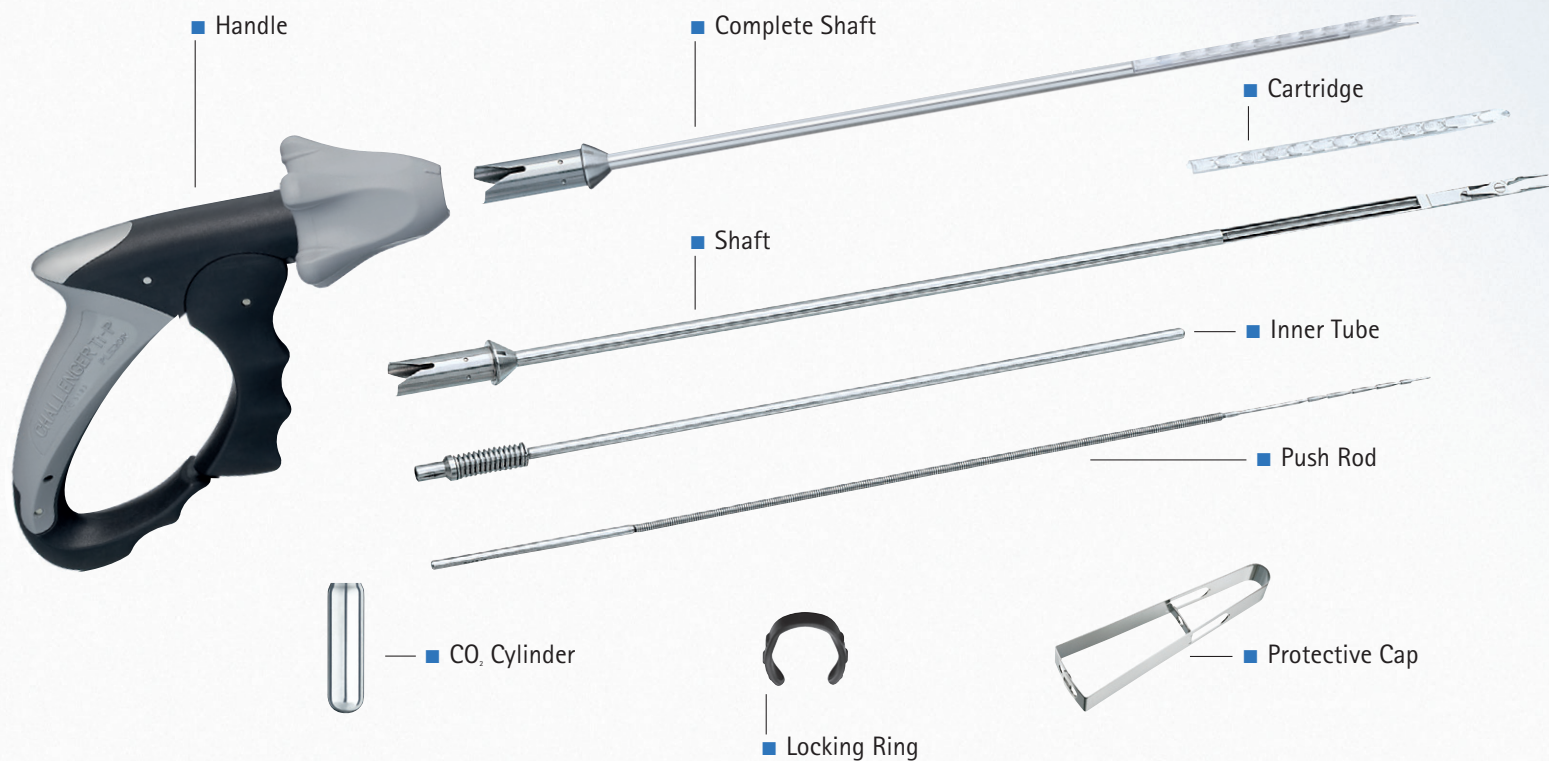

ML

SM

Spare Parts	Complete Shaft	PL536R	PL538R	PL522R
		Shaft Length 370 mm	Shaft Length 260 mm	Shaft Length 310 mm
		Diameter 10 mm	Diameter 10 mm	Diameter 5 mm
	Handle	PL520R	PL520R	PL520R
	Inner Tube	PL536802	PL538802	PL522801
		Shaft Length 370 mm	Shaft Length 260 mm	Shaft Length 310 mm
	Push Rod	PL536801	PL538801	PL522800
		Shaft Length 370 mm	Shaft Length 260 mm	Shaft Length 310 mm
	Locking Ring	PL520301	PL520301	PL520301
	CO ₂ Cylinder Only	PL575SU	PL575SU	PL575SU

Pack of 3, Single Sterile Packed

Preparation and Maintenance

O.R. Preparation

Ready to use in three steps!



Step 1
Assemble the handle and shaft



Step 2
Insert CO₂ cylinder




Step 3
Insert clip cartridge

Note: This product can be sterilized in an assembled or disassembled state. See *Instructions for Use* for more information.

Recommended products for care and handling

Protective Cap

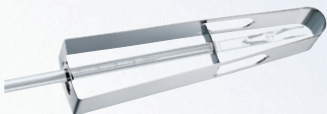

PL522810

Protects jaws during sterilization, washing and transportation
Diameter, 5 mm

PL536810

Protects jaws during sterilization, washing and transportation
Diameter, 10 mm

Sterilit® Oil Spray

JG600

For maintenance of instruments prior to sterilization
One unit

Cleaning Brushes



GK469R

Cleans inner tube diameter
5 mm, length 550 mm
Pack of 6 pcs.



PM995R

Cleans inner tube diameter
10 mm, length 550 mm
Pack of 6 pcs.

All rights reserved. Technical alterations are possible. This leaflet may be used for no other purposes than offering, buying and selling of our products. No part may be copied or reproduced in any form. In the case of misuse we retain the rights to recall our catalogs and price lists and to take legal actions.

©2016 AESCULAP. ALL RIGHTS RESERVED. PRINTED IN THE USA.
Aesculap is an equal opportunity employer

Aesculap, Inc. | 3773 Corporate Parkway | Center Valley, PA | 18034
Phone 800-282-9000 | Fax 610-791-6886 | www.aesculapusa.com

Aesculap, Inc. – a B. Braun company

DOC864 Rev A 11/16