

DS Ligation Clip System

Single-fire Clip Applier Solutions



Surgical Specialty Products

AESCULAP®

DS Ligation Clip System



Features & Benefits

Diamond-shaped clip texture


- Maximizes contact surface on tissue to further enhance clip security, while remaining atraumatic


"Double Clip" design

- Acts like a fenestrated grasper for a firm hold on tissue and to stabilize the clip against axial (side-to-side) dislodgement


Tooth & notch

- Stabilize the clip against radial (pull off) dislodgement

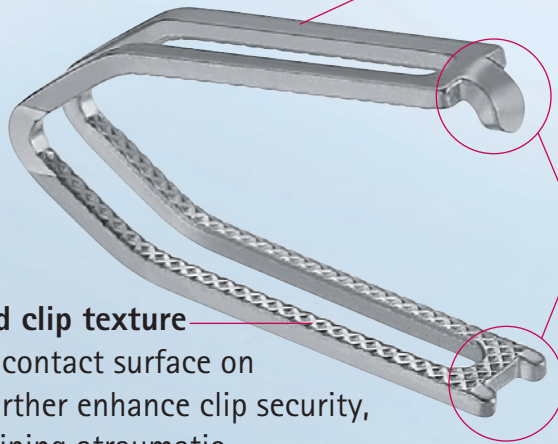

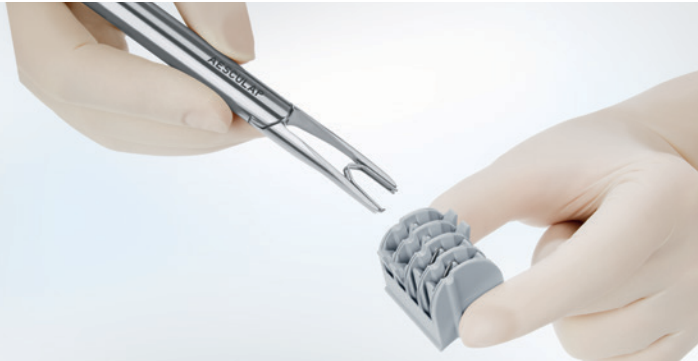






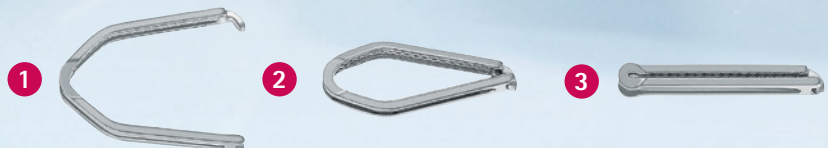



Closing behavior of the clip

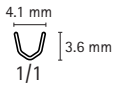


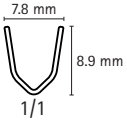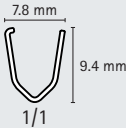





DS Clips

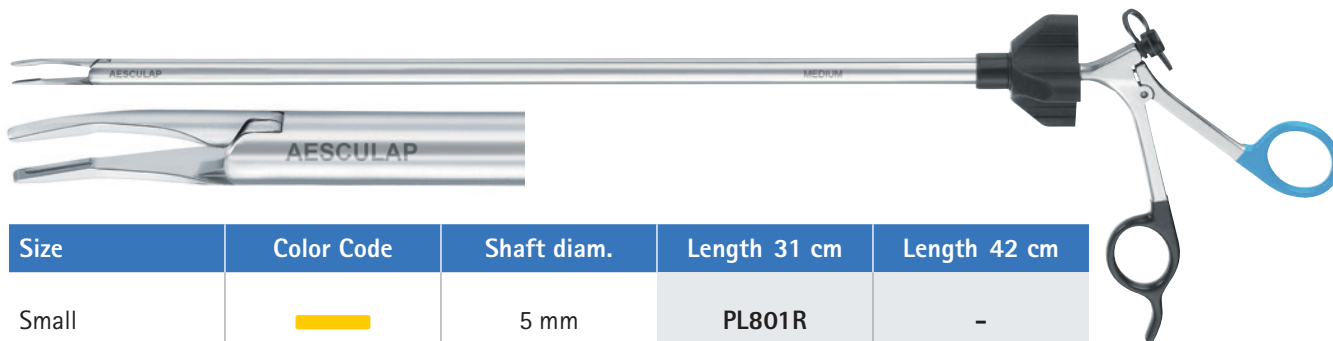Size	Color Code	Clip without Latch	Clip with Latch	Clips per Cartridge	No. of Cart. per Box
					
Small		 4.1 mm 3.6 mm 1/1	PL450SU	6	15
			PL452SU	24	6
Small-Medium		 5.0 mm 4.8 mm 1/1	PL453SU	6	15
Medium		 6.2 mm 5.8 mm 1/1	PL459SU	6	15
Medium-Large		 7.8 mm 8.9 mm 1/1	PL462SU	6	15
			 7.8 mm 9.4 mm 1/1		
Large		 10.3 mm 12 mm 1/1	PL468SU	6	12
			 10.3 mm 12.6 mm 1/1		
X-Large			 11.5 mm 16.2 mm 1/1	4	12


DS Clip Appliers and Removal Forceps




Laparoscopic Clip Appliers



Size	Color Code	Shaft diam.	Length 31 cm	Length 42 cm
Small		5 mm	PL801R	-
Small-Medium		5 mm	PL802R	PL822R
Medium		10 mm	PL806R	PL826R
Medium-Large		10 mm	PL807R	PL827R
Large		12 mm	PL808R	PL828R
X-Large		12 mm	PL809R	-

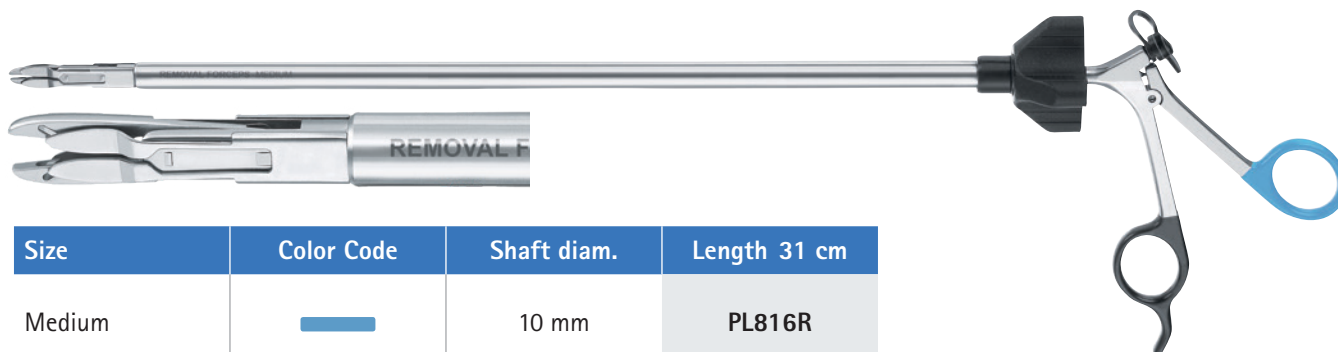Laparoscopic Clip Removal Forceps



Size	Color Code	Shaft diam.	Length 31 cm
Medium		10 mm	PL816R
Medium-Large		10 mm	PL817R
Large		12 mm	PL818R

Indication for Use:

Intended for use in endoscopic surgery for ligating and marking vessels and tubular structures whenever ligating clips are used/indicated.

All rights reserved. Technical alterations are possible. This leaflet may be used for no other purposes than offering, buying and selling of our products. No part may be copied or reproduced in any form. In the case of misuse we retain the rights to recall our catalogs and price lists and to take legal actions.

©2018 AESCULAP. ALL RIGHTS RESERVED. PRINTED IN THE USA.
Aesculap is an equal opportunity employer

Aesculap, Inc. | 3773 Corporate Parkway | Center Valley, PA | 18034
Phone 800-282-9000 | Fax 610-791-6886 | www.aesculapusa.com

Aesculap, Inc. – a B. Braun company

DOC1573 10/18