


Dual Mask

Dual Mask provides continual monitoring of carbon dioxide exhaled from both mouth and nose while delivering supplementary oxygen. Pre-fitted oxygen and sampling lines makes set-up quick and easy.

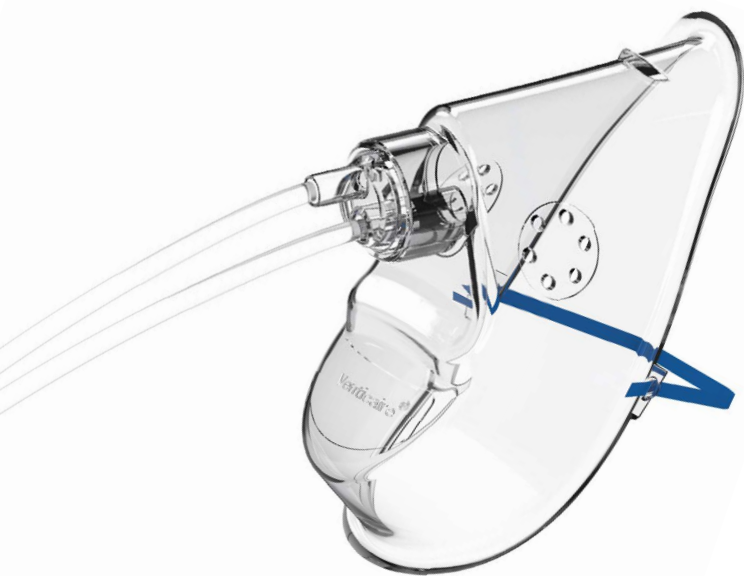

Easy Fit with Good Seal

The soft, pliable mask molds easily to the patient's face, with a feathered-edge rim providing maximum comfort.

Comfort and Security

Attached to the swivel connector, the oxygen tubing and monitoring line are lightly bonded together, but can be separated to reach the flow meter and capnometer, minimizing clutter around the patient.

Independently Tested

Dual Mask has been proven to be the most efficient in sampling end tidal CO₂ in an independent comparative study² with a simulated patient.

Star Lumen Tubing

Star lumen tubing reduces the risk of occlusion should the tube become kinked.

Excellent Waveforms

Unique design minimizes oxygen dilution, thus providing reliable sampling accuracy.



Male Luer



Female Luer



Microstream™ Luer

MALE LUER	FEMALE LUER	DESCRIPTION	TUBING	QTY
032-10-135U	032-10-134U	Adult Dual Mask with Nose Clip	10ft	50
032-12-143U	032-12-142U	Pediatric Dual Mask	10ft	50
032-10-137U	032-10-136U	Adult Dual Mask with Nose Clip and Filter	10ft	50
032-12-145U	032-12-144U	Pediatric Dual Mask and Filter	10ft	50

MICROSTREAM™ TYPE LUER	DESCRIPTION	TUBING	QTY
032-10-190U	Adult Dual Mask with Filter	7ft	50
032-10-191U	Adult Dual Mask with Filter	14ft	50
032-10-192U	Medium Dual Mask with Filter	7ft	50
032-12-190U	Pediatric Dual Mask with Filter	7ft	50
032-12-191U	Pediatric Dual Mask with Filter	14ft	50



Universal Adapter

LUER	DESCRIPTION	TUBING	QTY
032-10-135CU	Adult Dual Mask and Nose Clip with Male Luer + Universal Adapter	10ft	50
032-12-143CU	Pediatric Dual Mask and Nose Clip with Male Luer + Universal Adapter	10ft	50
032-10-141BU	Adult Dual Mask with Male Luer + Slip Disc Filter	10ft	50
032-12-145BU	Pediatric Dual Mask with Male Luer + Slip Disc Filter	10ft	50